
TiA 1.0

*Release***Christian Breitwieser, Christoph Eibel**

March 23, 2011

Contents

1	Signal Types	ii
1.1	Aperiodic Signals	ii
2	Control Connection	ii
2.1	Control Message Structure	ii
2.2	Reply Messages	iii
2.3	Server Commands	iii
	Check Protocol Version	iv
	Get MetaInfo	iv
	Get Data Transmission	iv
	Start Data Transmission	iv
	Stop Data Transmission	v
	Get Server State Connection	v
2.4	TiA Meta Info	v
2.5	TiA Error Description	v
3	Server State Connection	v
3.1	State Messages	v
	Server Running	vi
	Server Shut Down	vi
4	Data Packet	vii
4.1	Fixed Header	vii
	Byte Order	vii
	Packet Size	vii
	Time Stamp	vii
4.2	Variable Header	viii
4.3	Data	viii
5	Outlook	viii
5.1	Error Messages	viii
5.2	Datapackets	viii
6	Appendix	viii
6.1	TiA MetaInfo XML Schema	viii
6.2	TiA MetaInfo XML Example	x
6.3	TiA Error Message XML Schema	x

1 Signal Types

TiA defines the following signals types.

Identifier (string)	Value (hex)	Value (dec)	Description
eeg	0x00000001	1	Electroencephalogram
emg	0x00000002	2	Electromyogram
eog	0x00000004	4	Electrooculogram
ecg	0x00000008	8	Electrocardiogram
hr	0x00000010	16	Heart rate
bp	0x00000020	32	Blood pressure
button	0x00000040	64	Buttons (aperiodic)
joystick	0x00000080	128	Joystick axis (aperiodic)
sensors	0x00000100	256	Sensor
nirs	0x00000200	512	NIRS
fmri	0x00000400	1,024	FMRI
mouse	0x00000800	2,048	Mouse axis (aperiodic)
mouse-button	0x00001000	4,096	Mouse buttons (aperiodic)
not used yet			
user_1	0x00010000	65,536	User 1
user_2	0x00020000	131,072	User 2
user_3	0x00040000	262,144	User 3
user_4	0x00080000	524,288	User 4
undefined	0x00100000	1,048,576	undefined signal type
event	0x00200000	2,097,152	event

1.1 Aperiodic Signals

Aperiodic signals are signals which are not transmitted at constant time rates. For example the state of joystick axis.

Therefore values of the current state of these signals have be transmitted only if they have changed.

2 Control Connection

A TiA server has to provide a TCP port on which clients can create a control connection. Each client gets its own control connection to the server.

The control connection can be used by the clients to send requests to the server such as “start/stop data transmission”.

Some important remarks:

1. The protocol is in the style of HTTP (line structured text messages)
2. The messages are encoded in UTF-8
3. The message is split into lines which are terminated by the 0x0A character (also known as \n, “line feed” or <LF>)
4. Some messages contain additional XML-structured content which is UTF-8 encoded
5. All characters are case sensitive!

2.1 Control Message Structure

Each message which is send from the client to the server or vice versa is structured as followed:

1. Version line

2. Command line
3. Optional content description line
4. An empty line
5. Optional xml-structured content

Example:

```
TiA 1.0\n
CheckProtocolVersion\n
\n
```

Example with additional xml-structured content

```
TiA 1.0\n
MetaInfo\n
Content-Length: 79\n
\n
<?xml version="1.0" encoding="UTF-8"?><tiaMetaInfo version="1.0"></tiaMetaInfo>
```

2.2 Reply Messages

Each command message is answered with an reply message which contains either an OK or an Error. Error messages optionally contain an error description.

OK message:

```
TiA 1.0\n
OK\n
\n
```

Error message without an error description:

```
TiA 1.0\n
Error\n
\n
```

Error message including an error description:

```
TiA 1.0\n
Error\n
Content-Length: 73\n
\n
<tiaError version="1.0" description="Human readable error description."/>
```

2.3 Server Commands

A TiA 1.0 server implementation has to support the following commands:

- Check protocol version
- Get metainfo
- Get data connection
- Start data transmission
- Stop data transmission
- Get server state connection

Check Protocol Version

This command may be used by the client to check if the server understands the commands the client wants to send. The server has to respond with an OK message if it understands commands of the given protocol version.

Representation:

```
TiA 1.0\n
CheckProtocolVersion\n
\n
```

Server responses either with an OK or an error message.

Get MetaInfo

This command is used to get the informations about the signals from the server.

Representation:

```
TiA 1.0\n
GetMetaInfo\n
\n
```

Server response:

```
TiA 1.0\n
MetaInfo\n
Content-Length: [Length of XML Content in Bytes]\n
\n
<?xml version="1.0" encoding="UTF-8"?><tiaMetaInfo version="1.0">....</tiaMetaInfo>
```

or an error message.

Get Data Transmission

Two types of data transmissions exist: “TCP” and “UDP”.

Representation:

```
TiA 1.0 \n
GetDataConnection: TCP \n
\n
```

or

```
TiA 1.0 \n
GetDataConnection: UDP \n
\n
```

Server Response:

```
TiA 1.0 \n
DataConnectionPort: [Port-Number] \n
\n
```

or an error message.

Start Data Transmission

Representation:

```
TiA 1.0 \n
StartDataTransmission \n
\n
```

Server responses either with an OK or an error message.

Stop Data Transmission

Representation:

```
TiA 1.0 \n
StopDataTransmission \n
\n
```

Server responses either with an OK or an error message.

Get Server State Connection

Representation:

```
TiA 1.0\n
GetServerStateConnection\n
\n
```

Server Response:

```
TiA 1.0\n
ServerStateConnectionPort: [Port-Number]\n
\n
```

or an error message.

2.4 TiA Meta Info

The TiA meta info is structured in XML and contains information about the signals and the subject.

2.5 TiA Error Description

An error message in TiA version 1.0 optionally supports an error description in a human readable format. No error codes with special meaning are supported in this version.

3 Server State Connection

Additionally to the control- and data-connection a client may request a so called “server state connection”.

This connection is used by the server to transmit messages about its state to the clients.

3.1 State Messages

- Server running
- Server shut down

Server Running

Just to indicate that the server is running.

```
TiA 1.0\n
ServerStateRunning\n
\n
```

Clients must no reply to this message.

Server Shut Down


Indicates that the server will shut down soon. Therefore the control and the data connection will be closed by the server soon.

```
TiA 1.0\n
ServerStateShutdown\n
\n
```

Clients must not reply to this message.

4 Data Packet

TiA RAW Data Packet Structure:
(with exemplary content)

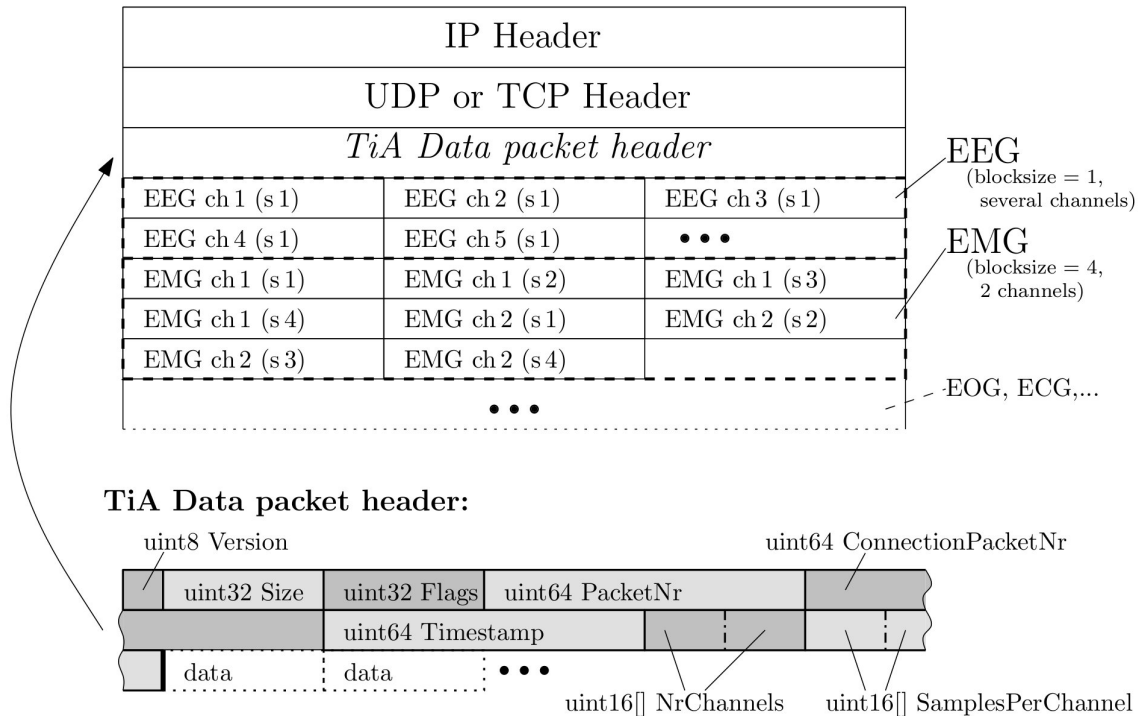

4.1 Fixed Header

Position (Byte)	Length (Byte)	Type	Content
0	1	unsigned integer	Datapacket Version has to be 3
1	4	unsigned integer	Packet size (in bytes), little endian
5	4	unsigned integer	Signal type flags
9	8	unsigned integer	Packet id
17	8	unsigned integer	Connection packet number
25	8	unsigned integer	Time stamp

Byte Order

All data and fields are transmitted in little endian.

Packet Size

The packet size is the total number of bytes (octets) of a datapacket including the fixed header, the variable header and the data.

The packet size of an empty datapacket is 33.

Time Stamp

microseconds since server start

4.2 Variable Header

Remarks:

- NoS is the number of signals. This number can be determined by counting the flags which are set in the “signal type flags” field of the fixed header.

Position (Byte)	Length (Byte)	Type	Content
33	2 * NoS	Array of 16bit unsigned integer	Number of channels
33 + (2 * NoS)	2 * NoS	Array of 16bit unsigned integer	Block size of signal (number of samples per channel)

4.3 Data

Data starts at the position $33 + (4 * \text{NoS})$. Sample values are transmitted as a sequence of 32bit floats (IEEE 754). Samples are firstly ordered by their signal type number (*Singal Types*).

5 Outlook

Future versions of TiA should guarantee downward compatibility. That means all TiA 1.0 commands should be supported by at least all 1.x versions.

5.1 Error Messages

Error messages may contain error codes to automatically interpret the meaning of an error.

5.2 Datapackets

Important remark: Due to downward compatibility no datapacket can become version 10! (as the first byte in datapacket version 2 is “10”).

6 Appendix

6.1 TiA MetaInfo XML Schema

```
1  <?xml version="1.0" encoding="UTF-8"?>
2  <!--
3      This is the XML Schema for the meta info XML representation version 1.0
4      of TOBI Interface A (TiA)
5
6      Version 1.0
7
8      author: Christoph Eibel
9  -->
10 <xsd:schema xmlns:xsd="http://www.w3.org/2001/XMLSchema">
11   <xsd:element name="tiaMetaInfo">
12     <xsd:complexType>
13       <xsd:sequence>
14
15         <!-- subject info -->
16         <xsd:element name="subject" minOccurs="0" maxOccurs="1">
```

```

17 <xsd:complexType>
18   <xsd:attribute name="id" type="xsd:string" use="optional"/>
19   <xsd:attribute name="firstName" type="xsd:string" use="optional"/>
20   <xsd:attribute name="surname" type="xsd:string" use="optional"/>
21   <xsd:attribute name="sex" use="optional">
22     <xsd:simpleType>
23       <xsd:restriction base="xsd:string">
24         <xsd:enumeration value="m" />
25         <xsd:enumeration value="f" />
26       </xsd:restriction>
27     </xsd:simpleType>
28   </xsd:attribute>
29   <xsd:attribute name="birthday" type="xsd:date" use="optional"/>
30   <xsd:attribute name="handedness" use="optional">
31     <xsd:simpleType>
32       <xsd:restriction base="xsd:string">
33         <xsd:enumeration value="l" />
34         <xsd:enumeration value="r" />
35       </xsd:restriction>
36     </xsd:simpleType>
37   </xsd:attribute>
38   <xsd:attribute name="medication" type="xsd:boolean" use="optional" />
39   <xsd:attribute name="glasses" type="xsd:boolean" use="optional" />
40   <xsd:attribute name="smoker" type="xsd:boolean" use="optional" />
41 </xsd:complexType>
42 </xsd:element>
43
44 <!-- master signal info -->
45 <xsd:element name="masterSignal" minOccurs="0" maxOccurs="1">
46   <xsd:complexType>
47     <xsd:attribute name="samplingRate" type="xsd:float" use="required"/>
48     <xsd:attribute name="blockSize" type="xsd:integer" use="required"/>
49   </xsd:complexType>
50 </xsd:element>
51
52
53 <!-- signal info -->
54 <xsd:element name="signal" minOccurs="0" maxOccurs="unbounded">
55   <xsd:complexType>
56
57     <!-- channels -->
58     <xsd:sequence>
59       <xsd:element name="channel" minOccurs="0" maxOccurs="unbounded">
60         <xsd:complexType>
61           <xsd:attribute name="nr" type="xsd:nonNegativeInteger" use="required"/>
62           <xsd:attribute name="label" type="xsd:string" use="required"/>
63         </xsd:complexType>
64       </xsd:element>
65     </xsd:sequence>
66
67     <!-- signal type: as defined in TiA 1.0 -->
68     <xsd:attribute name="type" use="required">
69       <xsd:simpleType>
70         <xsd:restriction base="xsd:string">
71           <xsd:enumeration value="eeg" />
72           <xsd:enumeration value="emg" />
73           <xsd:enumeration value="eog" />
74           <xsd:enumeration value="ecg" />
75           <xsd:enumeration value="hr" />
76           <xsd:enumeration value="bp" />
77           <xsd:enumeration value="buttons" />
78           <xsd:enumeration value="joystick" />
79           <xsd:enumeration value="sensors" />

```

```

80         <xsd:enumeration value="nirs" />
81         <xsd:enumeration value="fmri" />
82         <xsd:enumeration value="undefined" />
83     </xsd:restriction>
84 </xsd:simpleType>
85 </xsd:attribute>
86
87     <xsd:attribute name="samplingRate" type="xsd:float" use="required"/>
88     <xsd:attribute name="blockSize" type="xsd:integer" use="required"/>
89     <xsd:attribute name="numChannels" type="xsd:nonNegativeInteger" use="required"/>
90 </xsd:complexType>
91 </xsd:element>
92
93 </xsd:sequence>
94
95     <!-- version attribute -->
96     <xsd:attribute name="version" type="xsd:string" fixed="1.0" use="required" />
97
98 </xsd:complexType>
99 </xsd:element>
100
101 </xsd:schema>

```

6.2 TiA MetaInfo XML Example

```

1  <?xml version="1.0" encoding="UTF-8"?>
2  <tiaMetaInfo version="1.0">
3
4      <subject id="WE2" firstName="Max" lastName="Mustermann"
5          handedness="r" />
6
7      <masterSignal sampleRate="100" blockSize="10" />
8
9      <signal type="eeg" blockSize="10" sampleRate="100" numChannels="3">
10         <channel nr="1" label="Cz"/>
11         <channel nr="2" label="C1"/>
12         <channel nr="3" label="C2"/>
13     </signal>
14
15     <signal type="bp" blockSize="5" sampleRate="50" numChannels="5">
16         <channel nr="3" label="Channel 3 with Label"/>
17         <channel nr="2" label="Channel 2 with Label"/>
18     <!--
19         Channel 1, 4 and 5 have no labels
20     -->
21 </signal>
22
23 </tiaMetaInfo>

```

6.3 TiA Error Message XML Schema